



NAVIGATING THE SOLO₂[™] LCD SCREEN

QUICK REFERENCE GUIDE

NOTE: This guide is meant to be used only as a brief, general reference. For complete operating instructions, always refer to the Operator's Manual.

Turning Concentrator On/Off

NOTE: For this procedure, refer to FIGURE 1.1.

1. To turn the SOLO₂ concentrator on, press and hold the On/Off button  until the three warning lights display and the unit beeps.
2. To turn the SOLO₂ concentrator off, press and hold the On/Off button .

NOTE: GREEN light will flash and you will hear the machine stop running.

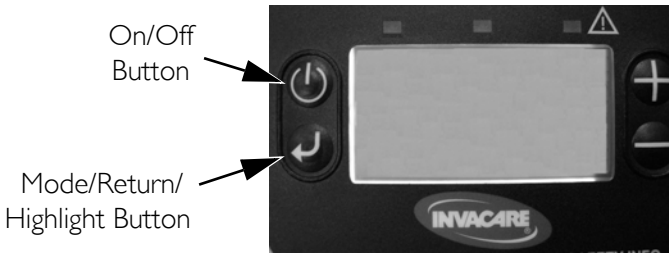



FIGURE 1.1 Turning Concentrator On/Off

Menu Screens Navigation

NOTE: For this procedure, refer to FIGURE 1.1.

1. After turning on unit, press and hold the Mode/Return/Highlight Button  to scroll through the three SOLO₂ concentrator menu screens.

NOTE: Screens will continue scrolling until the Mode/Return/Highlight Button is released.

Normal Operating/Standby Screen

The Normal Operating/Standby screen shows either battery strength gauge, fan operating, external power on and battery charging icons, when applicable, or the current/last mode of operation and output flow level selected.

LCD Adjust Screen







The LCD Adjust screen allows user to change the brightness and contrast of the screen (FIGURE 1.2).



FIGURE 1.2 LCD Adjust Screen

LCD Adjust Screen Navigation

NOTE: For this procedure, refer to FIGURE 1.2.

1. Perform one of the following, while the LCD Adjust screen is displayed:
 - A. Tap the Mode/Return/Highlight Button  to move the black highlight bar from left to right, or right to left.
 - B. Tap the Up/Increase Button  or the Down/Decrease Button  while the black highlight bar is on the “LEFT SIDE” of the screen to move the highlight bar up or down.
 - C. Tap the Up/Increase Button  or Down/Decrease Button  while the highlight bar is on the “RIGHT SIDE” of the screen to change the value of the item highlighted.
 - D. Hold down the Mode/Return/Highlight Button  to continue scrolling through the three menu screens.



NOTE: If no buttons are pressed for approximately 10 seconds, the screen will change back to the Standby/Normal Operating screen. Any changes made will NOT be saved.

2. To save changes made on the LCD Adjust screen:

- A. Press and hold the Mode/Return/Highlight Button  until a different screen appears.

Mode & Flow Adjust Screen

The Mode & Flow Adjust screen allows user to change the current mode of operation and output flow (FIGURE 1.3).

NOTE: To adjust either the mode or flow, the Mode & Flow Adjust screen must be displayed. If the battery level screen is displayed, press the Mode/Return/Highlight , Up/Increase , or Down/Decrease button for approximately two seconds (screen should change to Mode & Flow Adjust).


NOTE: Any changes to the setting will not take effect until after the screen flashes three times.



FIGURE 1.3 LCD Adjust Screen



Changing Mode

NOTE: For this procedure, refer to FIGURE 1.3.

1. With the Mode & Flow screen displayed, press and hold the Mode/Return/Highlight Button  for one second.
2. Wait until the Mode & Flow screen flashes three times.

Adjusting Flow

NOTE: For this procedure, refer to FIGURE 1.3.

1. With the Mode & Flow screen displayed, press and hold the Up/Increase Button  or Down/Decrease Button  until desired setting is achieved.
2. Wait until the Mode & Flow screen flashes three times.



Yes, you can.

Invacare Corporation

www.invacare.com

USA

One Invacare Way
Elyria, Ohio USA
44036-2125
800-333-6900
Technical Services
800-832-4707

© 2010 Invacare Corporation. All rights reserved. Republication, duplication or modification in whole or in part is prohibited without prior written permission from Invacare. Trademarks are identified by [™] and [®]. All trademarks are owned by or licensed to Invacare Corporation or its subsidiaries unless otherwise noted.

Part No 1163220

Rev B - 01/10

